

Balloon Trocars

For atraumatic and secure cannula fixation during laparoscopic surgery

Balloon Trocars

The latex-free balloon trocars from LaproSurge provide enhanced stability throughout laparoscopic procedures by additional fixation from above and below the abdominal wall, reducing unintentional cannula migration.

The atraumatic stabilisation helps to minimise fascial trauma while also maintaining pneumoperitoneum, even after incisions have been enlarged for specimen removal.

After trocar insertion, the balloon is inflated with air through the cannula's inflation port using the syringe supplied.

Once in position, the gel anchor or retention disk is adjusted to the required cannula depth to gently fix the cannula to the patient's abdomen.




- 5mm, 10/11mm and 12mm balloon trocars supplied with a retention disk for stability
- 12mm balloon trocar also available with an adjustable gel anchor grip with suture tabs for additional fixation




Reduced Abdominal Intrusion

5mm and 10mm-12mm balloon trocars reduce abdominal intrusion by 48% and 44% respectively compared to standard threaded cannulas by fixing the cannula into position at the distal end.

The reduced intra-abdominal cannula depth increases working space and in turn enhances visibility.

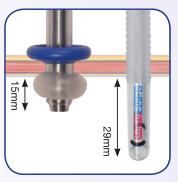

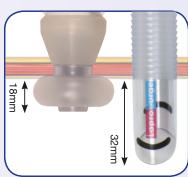

Ordering Code	Fixation	Length	Box Qty
BC5SS-C	Retention disk	70mm	5 units
BC5SS-O	Retention disk	70mm	5 units
BC5SL-C	Retention disk	100mm	5 units
BC5SL-O	Retention disk	100mm	5 units
BC5SL	Retention disk	100mm	5 units
BC10SLR-C Auto	Retention disk	100mm	5 units
BC10SLR-O Auto	Retention disk	100mm	5 units
BC10SLR Auto	Retention disk	100mm	5 units
BC12SL-C Auto	Gel cone	100mm	5 units
BC12SHL Auto	Gel cone	100mm	5 units
BC12SL-O Auto	Gel cone	100mm	5 units
BC12SLR-C Auto	Retention disk	100mm	5 units
BC12SLR-O Auto	Retention disk	100mm	5 units
	BC5SS-C BC5SS-O BC5SL-C BC5SL-O BC5SL BC10SLR-C Auto BC10SLR-O Auto BC10SLR Auto BC12SL-C Auto BC12SHL Auto BC12SHL Auto BC12SL-O Auto	BC5SS-C Retention disk BC5SS-O Retention disk BC5SL-C Retention disk BC5SL-O Retention disk BC5SL Retention disk BC10SLR-C Auto Retention disk BC10SLR-O Auto Retention disk BC10SLR Auto Retention disk BC12SL-C Auto Gel cone BC12SHL Auto Gel cone BC12SL-O Auto Gel cone BC12SLR-C Auto Retention disk	BC5SS-O Retention disk 70mm BC5SL-C Retention disk 100mm BC5SL Retention disk 100mm BC10SLR-C Auto Retention disk 100mm BC10SLR Auto Retention disk 100mm BC12SL-C Auto Retention disk 100mm BC12SHL Auto Gel cone 100mm BC12SL-O Auto Gel cone 100mm BC12SL-O Auto Gel cone 100mm BC12SL-O Auto Retention disk 100mm



For more information visit our website or contact us:

Tel: +44 (0)1923 220 711 Fax: +44 (0)20 8950 8659 E-mail: info@laprosurge.com