




Lapro Surge®

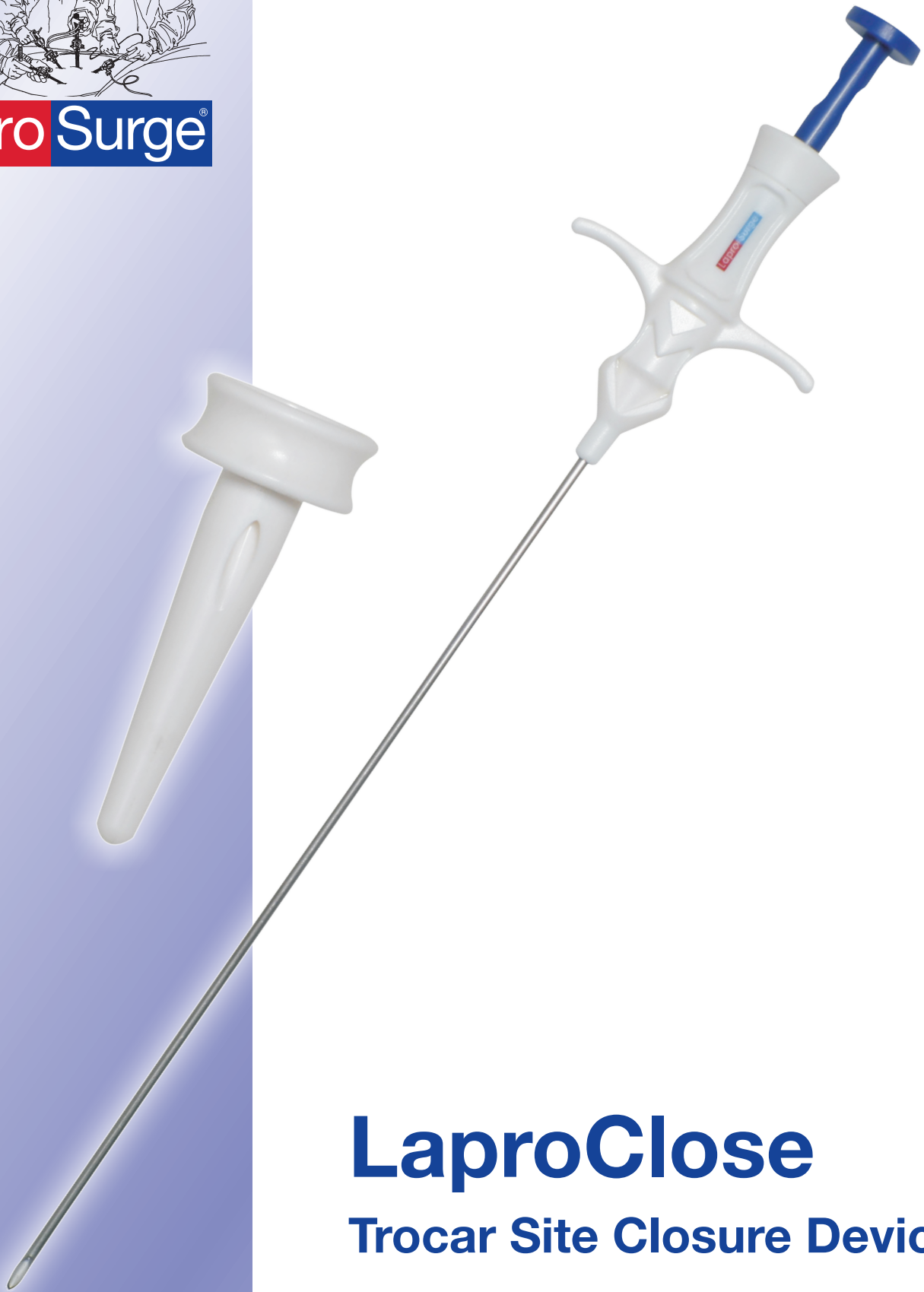

# LaproClose

## Trocar Site Closure Device

---

A simple and reliable system for fascial closure of port site defects

# LaproClose Trocar Site Closure Device

Complete system for quick, safe and reliable closure of laparoscopic incision sites.

## Cone-shaped guide for closure of 10mm/12mm incisions


- Correctly angles the introducer to capture all tissue layers for full closure of fascia and peritoneum
- Prevents loss of pneumoperitoneum during port site closure



## Rigid introducer for passing sutures through all tissue layers

Available with or without guide.

Intuitive spring mechanism for extending and retracting the grasping jaws

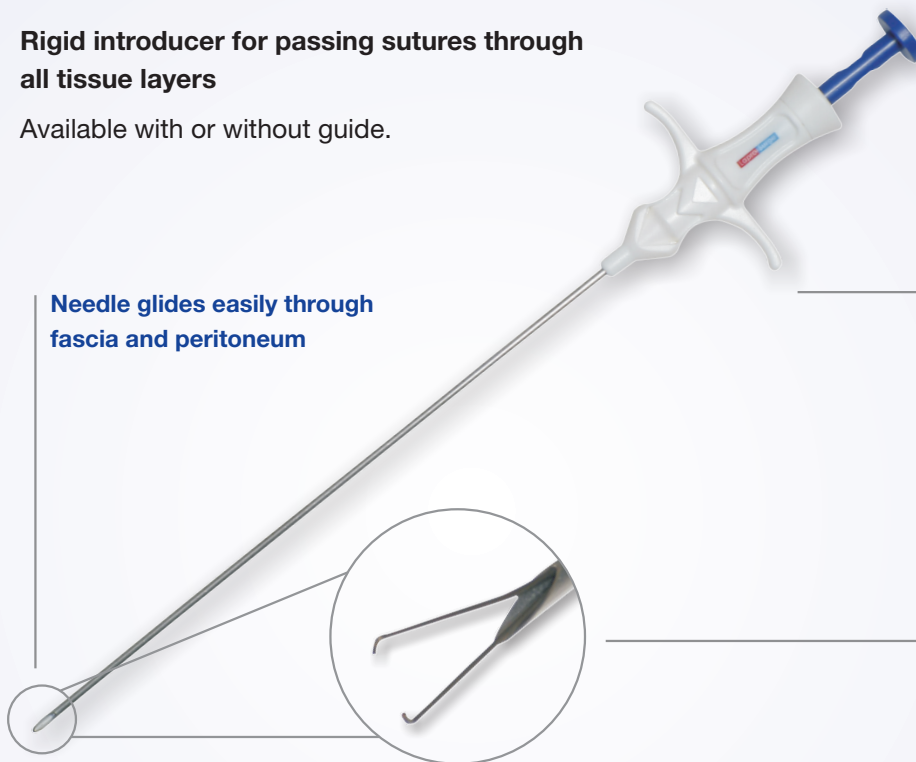

Needle glides easily through fascia and peritoneum

Ergonomic handle for comfortable finger placement

Long, retractable grasping jaws to accommodate any suture size

The combined use of the guide and introducer together minimises the risk of internal injuries by ensuring the introducer passes through the abdominal wall in an acute angle, enabling better visualisation of the needle entering the peritoneal cavity.

Product	Ordering Code	Box Qty
Complete Device – Guide and Introducer	TSC001	10
Introducer only	TSC002	10



For more information visit our website or contact us:

Tel: +44 (0)20 8950 8662 Fax: +44 (0)20 8950 8659 E-mail: [info@laprosurge.com](mailto:info@laprosurge.com)